

3.2.1 BeamGage

See your beam as never before

3.2.1.1 BeamGage®-Standard Version

- Extensive set of ISO quantitative measurements
- Patented Ultracal™ algorithm for highest accuracy measurements in the industry
- Customizable user interface for 'ease of use'
- Auto-setup and Auto-exposure capabilities for fast set-up and optimized accuracy
- Statistical analysis on all calculated results displayed in real time
- New BeamMaker® beam simulator for algorithm self-validation

The performance of today's laser systems can strongly affect the success of demanding, modern laser applications.

The beam's size, shape, uniformity or approximation to the expected power distribution, as well as its divergence and mode content can make or break an application. Accurate knowledge of these parameters is essential to the success of any laser-based endeavor. As laser applications push the boundaries of laser performance it is becoming more critical to understand the operating criteria.


For over thirty years Ophir-Spiricon has developed instruments to accurately measure critical laser parameters. Our LBA and BeamStar software have led the way. Now with the introduction of BeamGage, Ophir-Spiricon offers the first "new from the ground up" beam profile analysis instrument the industry has experienced in over 10 years.

BeamGage includes all of the accuracy and ISO approved quantitative results that made our LBA software so successful. BeamGage also brings the ease-of-use that has made our BeamStar software so popular. Our patented UltraCal algorithm, guarantees the data baseline or "zero-reference point" is accurate to 1/10 of a digital count on a pixel-by-pixel basis. ISO 11146 requires that a baseline correction algorithm be used to improve the accuracy of beam width measurements. UltraCal has been enhanced in BeamGage to assure that accurate spatial measurements are now more quickly available.




See Your Beam As Never Before:


The Graphical User Interface (GUI) of BeamGage is new. Dockable and floatable windows plus concealable ribbon tool bars empowers the BeamGage user to make the most of a small laptop display or a large, multi-monitor desktop PC.

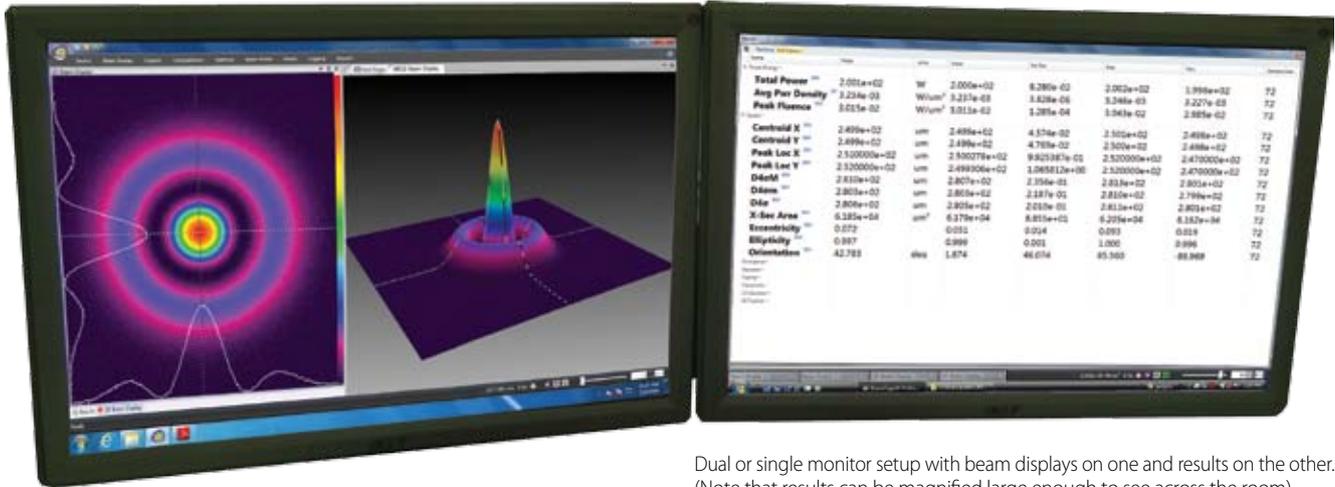


Dual or single monitor setup with beam displays on one and results on the other. (Note that results can be magnified large enough to see across the room).

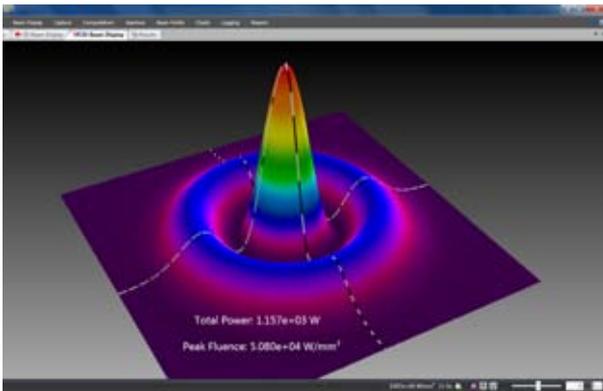

Beam only (Note results overlaid on beam profile).

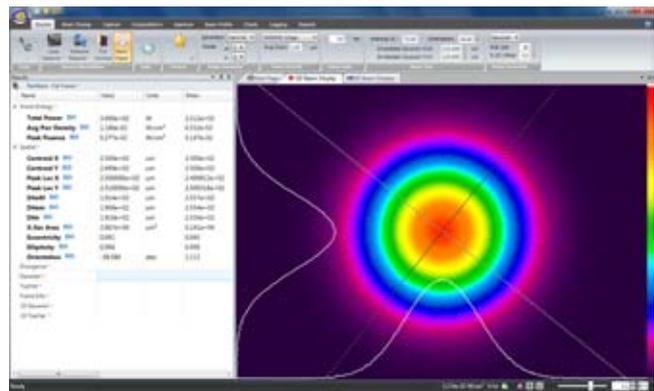

Beam plus results

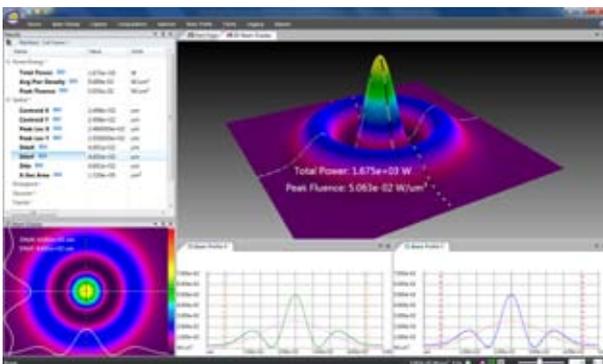

Multiple beam and results windows. (Note quantified profile results on 3D display & quantified 2D slices).

- 3D displays Rotate & Tilt. All displays Pan, Zoom, Translate & Z axis Zoom.

Measure Your Beam As Never Before:

Ultracal: Essential, or no big deal?


If you want accurate beam measurements, you want Ultracal.

What is Ultracal?


Our patented, baseline correction algorithm helped establish the ISO 11146-3 standard for beam measurement accuracy. The problems with cameras used in beam profile measurements are: a) baseline, or zero, of the cameras drift with time temperature, and b) include random noise. Ultracal is the only beam profiler algorithm that sets the baseline to “zero”, and, in the center of the noise. (Competitive products use other less sophisticated algorithms that perform a baseline subtraction, but truncate the noise below the “zero” of the baseline. This leaves only a “positive” component, which adds a net value to all beam measurements).

Try the following on any other beam profiler product to see the inherent error if you don't use Ultracal.

1. Measure a beam with full intensity on the profiler camera.
2. Insert a ND2 filter (100X attenuation) into the beam and measure it again.
3. Compare the results.
4. The Standard Deviation below is about 3%, which is phenomenal compared to the 100% or more of any beam profiler without Ultracal.

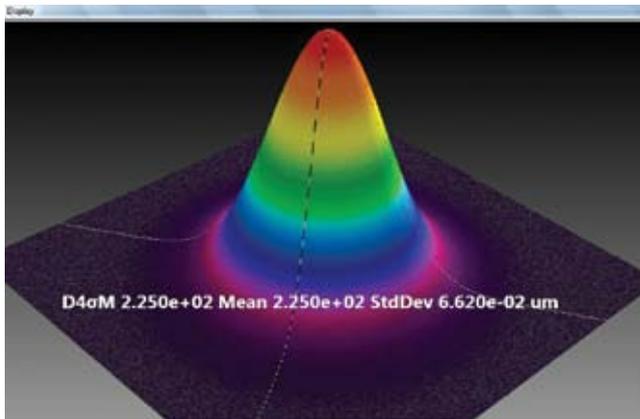

Beam at full intensity, Width 225 μ m, Std Dev 0.06 μ m

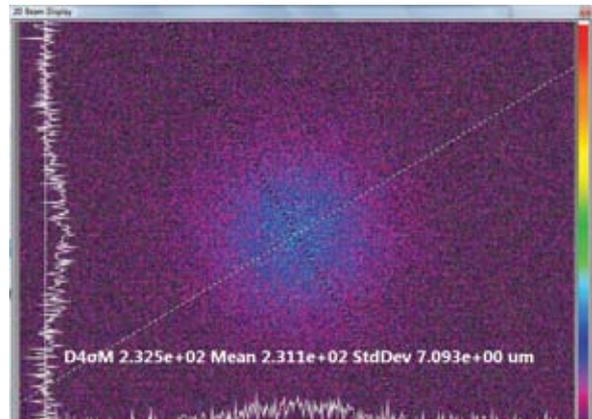

Beam attenuated 100X (displayed here in 2D at 16X magnitude zoom), Width 231 μ m, Std Dev 7 μ m

Adding the use of Automatic Aperture improves the accuracy to 1%. (The conditions of this measurement is a camera with a 50dB SNR).

5. You normally don't make measurements at such a low intensity. But occasionally you may have a drop in intensity of your beam and don't want to have to adjust the attenuation. Or, you may occasionally have a very small beam of only a few tens of pixels. In both of these cases, Ultracal becomes essential in obtaining accurate measurements.

Beam Measurements and Statistics

BeamGage allows you to configure as many measurements as needed to support your work, and comes standard with over 55 separate measurement choices. To distinguish between calculations that are based on ISO standards and those that are not, a graphical ISO logo is displayed next to appropriate measurements. You can also choose to perform statistical calculations on any parameter in the list.

Name	Value	Units
Spatial *		
Centroid X ISO	3.121e+00	mm
Centroid Y ISO	3.121e+00	mm
Peak Loc X ISO	3.10000e+00	mm
Peak Loc Y ISO	3.12500e+00	mm
D4σM ISO	4.449e+00	mm
D4σm ISO	4.406e+00	mm
Dkcm 10/90	3.778e+00	mm
Dkem 10/90	3.685e+00	mm
Dkcm 16/84 ISO	3.477e+00	mm
Dkem 16/84 ISO	3.368e+00	mm
D%pkM	2.714e+00	mm
D%pkm	2.594e+00	mm
X-Sec Area ISO	1.540e+01	mm ²
Eccentricity ISO	0.138	
Ellipticity ISO	0.990	
Divergence *		
Gaussian *		
Gauss Centroid X	3.125039e+00	mm
Gauss Centroid Y	3.124977e+00	mm
Goodness of Fit	0.694	
Roughness of Fit	0.217	


Small sample of possible measurements out of a list of 55

Name	Value	Units	Mean	Std Dev	Max	Min	Sample Size
Power/Energy *							
Total Power ISO	2.809e+02	W	2.809e+02	4.096e-02	2.810e+02	2.808e+02	248
Peak Fluence ISO	8.105e+01	W/mr	8.111e+01	1.558e-01	8.170e+01	8.073e+01	248
Efficiency ISO	---	%	---	---	---	---	---
% in Aperture	100.00	%	100.00	0.00	100.00	100.00	248
Spatial *							
Centroid X ISO	3.122e+00	mm	3.121e+00	2.820e-04	3.122e+00	3.121e+00	262
Centroid Y ISO	3.122e+00	mm	3.121e+00	2.683e-04	3.122e+00	3.121e+00	262
Peak Loc X ISO	3.125000e+00	mm	3.124046e+00	2.578e-02	3.200000e+00	3.050000e+00	262
Peak Loc Y ISO	3.125000e+00	mm	3.128721e+00	2.567e-02	3.200000e+00	3.075000e+00	262
D4σM ISO	4.451e+00	mm	4.450e+00	1.176e-03	4.454e+00	4.435e+00	1.733
D4σm ISO	4.406e+00	mm	4.407e+00	1.208e-03	4.421e+00	4.403e+00	1.733
Dkcm 10/90	3.767e+00	mm	3.770e+00	5.985e-03	3.788e+00	3.750e+00	262
Dkem 10/90	3.674e+00	mm	3.676e+00	6.629e-03	3.695e+00	3.653e+00	262
Eccentricity ISO	0.141		0.139	0.003	0.147	0.132	262
Ellipticity ISO	0.990		0.990	0.000	0.991	0.989	48
Divergence *							
Gaussian *							
TopHat *							
Frame Info *							
1D Gaussian *							
1D TopHat *							

Sample of calculation results with statistics applied

Multiple Charting Options

You can create strip charts for stability observations on practically any of the calculations options available. Charts enable tracking of short or long term stability of your laser.

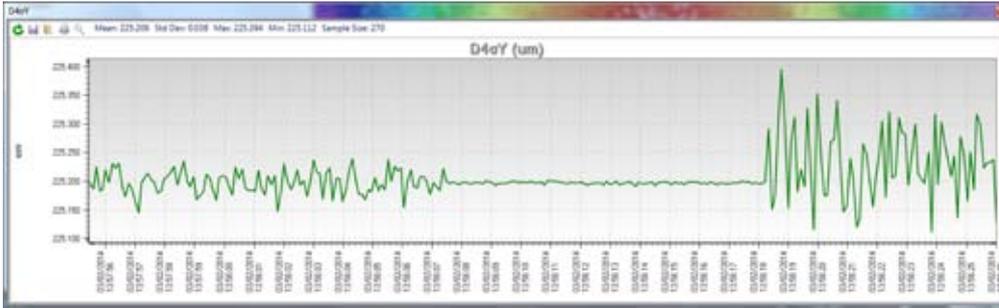

Strip chart of beam D4sigma width. Note how changing conditions affects the width repeatability. Beam intensity changed over 10db, making noise a significant factor in measurement stability.

Beam Pointing Stability

Open the Pointing Stability Window to collect centroid and peak data from the core system and display it graphically. View a chart recorder and statistical functions in one interface:

The screenshot shows the BeamGage software interface with several windows open. On the left, a 'Peak Energy' table lists statistics like Total, Peak, and Min. In the center, a 'Centroid X (um)' strip chart shows data over time. On the right, a 'Centroid Y (um)' strip chart shows similar data. A 'Peak location scatter plot with histogram color-coding' is visible in the top right, and a 'Centroid location scatter plot with histogram color-coding' is in the middle right. Red callout boxes point to these features with descriptive text.

Peak location scatter plot with histogram color-coding.

Set a sample limit, and specify the results items to graph on the strip chart.

The radius is referenced from either an Origin established in BeamGage or from the continuously calculated Average Centroid position.

A centroid location scatter plot with histogram color-coding.

A pointing stability strip chart presents data over time for the Centroid X and Y, Peak X and Y and centroid radius from an origin or from the mean centroid.

Easy to Use and Powerful

BeamGage is the only beam profiler on the market using modern Windows 7 navigation tools. The menu system of BeamGage is easy to learn and easy to use with most controls only one mouse click away. Some ribbon toolbar examples:



Some of the Beam Display options (Display access options under the Tools tab on the left).



Some of the Beam Capture options.

BeamGage Main Display Screen

The screenshot shows the BeamGage software interface with several callout boxes pointing to specific features:

- File Save/Load ApplicationButton**: Points to the top-left corner of the application window.
- Quick Access Toolbar for common tasks**: Points to the toolbar below the menu bar.
- Tabbed Control Access**: Points to the tabs for 'Computations', 'Aperture', 'Beam Profile', 'Charts', 'Logging', and 'Reports'.
- 2D Beam Display**: Points to the central 2D beam profile plot.
- Tool Windows that dock inside or float outside App**: Points to the 'Results' and 'D4σ' tool windows.
- User Definable Window Layout**: Points to the overall arrangement of the software windows.
- Integrated Help System**: Points to the help icon in the top-right corner.
- Beam Results With Statistics**: Points to the 'Results' table on the left.
- ISO Compliant Results**: Points to the 'ISO' column in the results table.
- 1D Profiling Options**: Points to the 'Beam Profile X' and 'Beam Profile Y' plots.
- Cursors With Power / Energy Readouts**: Points to the vertical cursors on the 1D profiles.
- Processing Status Indicators**: Points to the status bar at the bottom.
- 3D Beam Display**: Points to the 3D beam profile plot at the bottom right.
- Buffered Video Scrolling Controls**: Points to the video playback controls at the bottom right.

Pass / Fail with Password Protection for Production Testing

BeamGage allows the user to configure the displayed calculations; set-up the screen layout and password protect the configuration from any changes. This permits secure product testing as well as data collection for Statistical Process Control (SPC), all while assuring the validity of the data.

Name	Value	Mean	Std Dev	Max	Min	Units
Power/Energy *						
Spatial *						
Centroid X ISO	7.831e+01	7.831e+01	2.849e-03	7.832e+01	7.830e+01	um
Centroid Y ISO	7.965e+01	7.965e+01	3.047e-03	7.966e+01	7.964e+01	um
Peak Loc X ISO	7.000000e+01	7.073199e+01	1.340173e+0	7.500000e+01	6.700000e+01	um
Peak Loc Y ISO	7.100000e+01	7.183659e+01	1.333245e+0	7.500000e+01	6.800000e+01	um
D4σX ISO	1.238e+02	1.238e+02	8.334e-03	1.239e+02	1.238e+02	um
D4σY ISO	124.041	124.053	0.008	124.079	124.027	um
D4σ ISO	1.239e+02	1.239e+02	6.395e-03	1.240e+02	1.239e+02	um

Failures (or successes) can be the impetus for additional actions including a TTL output signal or PC beep and the termination of further data acquisition.

Camera Compatibility

For lasers between 190-1100nm wavelengths, BeamGage interfaces to silicon CCD USB and GigE cameras. For applications between 1440-1605nm, BeamGage supports cost effective phosphor coated CCD cameras. For demanding applications between 900-1700nm, BeamGage supports an InGaAs camera. And for applications in the ultraviolet, 13-355nm, or far infrared or Terahertz range, 1.06-3000nm, BeamGage supports Spiricon's Pyrocam, pyroelectric array cameras.

190-1100nm*



Model	SP503U	SP620U	SP300	Gevicam
Spectral Response nm	190 - 1100nm*	190 - 1100nm*	190 - 1100nm*	190 - 1100nm*
Application	1/2" format, slim profile, wide dynamic range, CW & pulsed lasers, adjustable ROI	1/1.8" format, high resolution, wide dynamic range, CW & pulsed lasers, adjustable ROI	1/1.8" format, high resolution, CW & pulsed lasers, adjustable ROI	1/1.8" format, high resolution, networkable, long cable distances, adjustable ROI
Number of Elements	640 x 480	1600 x 1200	1928 x 1448	1600 x 1200
Interface Style	USB 2.0	USB 2.0	USB 3.0, USB 2.0	GigE
Windows OS support	Windows 7 (32/64)			

190-1100nm*



Model	LT665	L11059
Spectral Response nm	190 - 1100nm*	190 - 1100nm*
Application	12.5mm x 10mm, 1" format for large beams, CW & pulsed lasers, adjustable ROI	36mm x 24mm, 35mm format for large beams, CW & pulsed lasers, adjustable ROI
Number of Elements	2752 x 2192	4008 x 2672
Interface Style	USB 3.0, USB 2.0	USB 2.0
Windows OS support	Windows 7 (32/64)	

*Although our silicon cameras have shown response out to 1320nm it can cause significant blooming which could lead to significant errors of beam width measurements. We would suggest our XC13 InGaAs camera for these wavelengths to give you the best measurements.

1440-1605nm



Model	SP503U-1550	SP620U-1550	LT665-1550
Spectral Response nm	1440 - 1605nm	1440 - 1605nm	1440 - 1605nm
Application	NIR wavelengths, 1/2" format, low resolution	NIR wavelengths, 1/1.8" format, low resolution, adjustable ROI and binning	12.5mm x 10mm, 1" format for large beams, CW & pulsed lasers, adjustable ROI
Number of Elements	640 x 480	1600 x 1200	2752 x 2195
Interface Style	USB 2.0	USB 2.0	USB 3.0, USB 2.0
Windows OS support	Windows 7 (32/64)		

900-1700nm



Model	XEVA 100Hz
Spectral Response nm	900 - 1700nm
Application	High resolution InGaAS performance, NIR wavelengths
Number of Elements	320 x 256
Interface Style	USB 2.0
Windows OS support	Windows 7 (32)

13-355nm & 1.06-3000µm



Model	Pyrocam IIIHR	Pyrocam IV
Spectral Response nm	13-355nm & 1.06-3000µm	13-355nm & 1.06-3000µm
Application	UV & Far IR Only commercial array to view Terahertz	UV & Far IR Only commercial array to view Terahertz
Number of Elements	160 x 160	320 x 320
Interface Style	GigE	GigE
Windows OS support		Windows 7 (32/64)

Unique Features of BeamGage - Standard

Power/Energy Calibration


Using the USB output from select Ophir power/energy meters, the BeamGage application will display measured power/energy values from the full range of Ophir thermopile, photodiode and pyroelectric sensors. Pulsed lasers can be synced up to 100Hz, or the frame rate of the triggered camera, whichever is less. This is the first time in the industry a laser power meter has been married to a laser beam profile system.



BeamGage is the only product to integrate profiling and power meter measurements

BeamMaker®; Numerical Beam Profile Generator

BeamGage contains a utility, BeamMaker, that can synthetically generate beam profile data by modeling either Laguerre, Hermite or donut laser beams in various modal configurations. BeamMaker permits the user to model a beam profile by specifying the mode, size, width, height, intensity, angle, and noise content. Once generated the user can then compare the theoretically derived measurements to measurements including experimental inaccuracies produced by the various measurement instruments and environmental test conditions. Users can now analyze expected results and confirm if measurement algorithms will accurately measure the beam even before the experiment is constructed. BeamMaker can help laser engineers, technicians and researchers understand a beam's modal content by calculating results on modeled beams for a better understanding of real laser beam profiles. BeamMaker is to laser beam analysis as a function generator is to an oscilloscope.

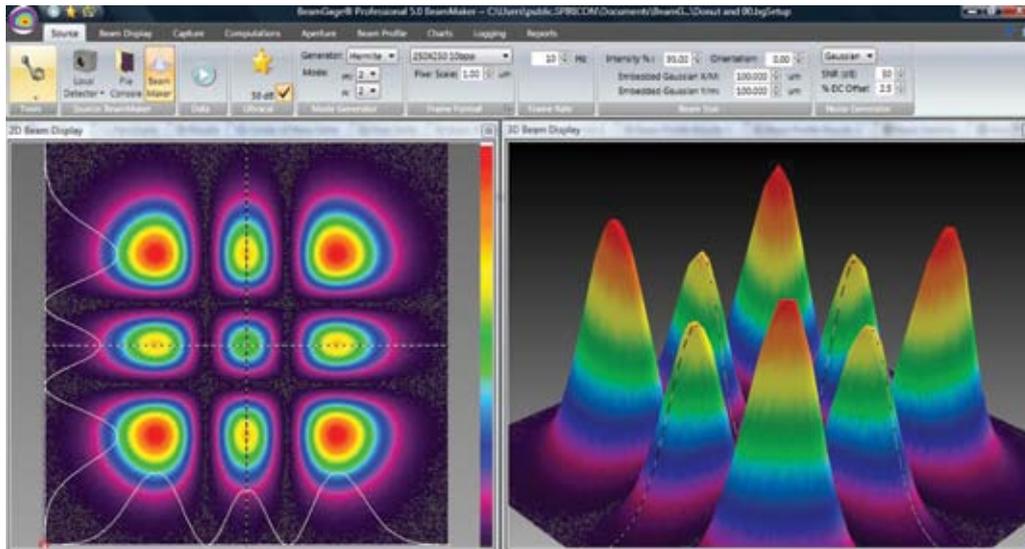

BeamMaker producing a synthetically generated Hermite TEM₂₂ beam and displayed in both 2D and 3D

Integrated automatic Help linked into the Users Guide

Touch sensitive Tool tips are available on most all controls, and "What's This" help can provide additional details. Confused about what something is or forgot how it works, just go to the top right corner and touch the "What's This" help icon, then click on the control or menu item that you want more info about and you are taken to the explanation within the BeamGage Users Guide.

Multilingual


BeamGage comes with both Japanese and Chinese user interface. Country specific manuals can be downloaded from the ophiropt.com/photronics web site.

3.2.1.2 BeamGage®-Professional Version

Professional is an upgrade version of BeamGage-Standard that has all of the BeamGage-Standard features plus additional functionality.

Image Partitioning

Partitioning allows the user to subdivide the camera image into separate regions, called partitions, and compute separate beam results within each partition. When using partitioning special results items can be displayed that relate to delta values between the computed centroids or peaks of each partition. Partitioning is useful to enable separate analysis of individual beams when multiple beams impinge on the camera simultaneously. This feature is particularly useful when analyzing multiple fibers in a single bundle.

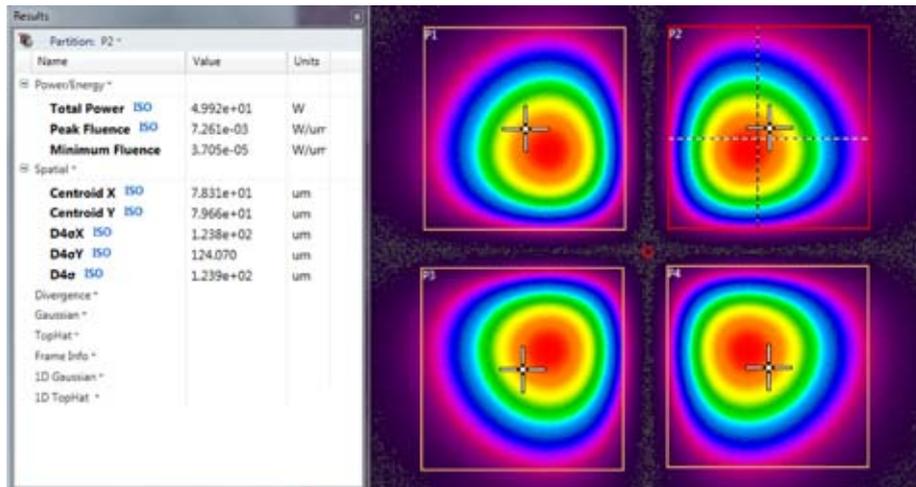

Shown is an example of the results for partition P2 and its related display frame. Observe that the selected partition is highlighted in RED. The crosshair in each partition is user controlled. The crosshair can be moved to a new position with the mouse or can be numerically positioned using the expanded controls that appear when a partition is created.

Automation Interface

BeamGage Professional provides an automation interface via .NET components to allow customers the ability to build custom applications that incorporate the laser beam analysis and processing power of BeamGage. The BeamGage automation interface allows developers to control BeamGage programmatically via a set of "puppet strings" known as the automation interface. The automation interface was developed to provide the ability to base control decisions for a second application on results and behaviors recognized by BeamGage. With this ability users can quickly and efficiently meet their manufacturing/analysis goals with minimum human interaction.

The automation interface was designed to achieve two main goals. First, to allow the BeamGage user to programmatically do what they could otherwise do via the graphical user interface (GUI). Second, to expose stable interfaces to the user that will not change, causing breaks to their dependant code. Interface examples for LabVIEW, Excel and .NET VB are included.

Custom Calculations


If BeamGage-Standard does not have the measurement you need the Professional and Enterprise versions permit the user to program-in their own set of calculations. User defined computations are treated the same as other BeamGage standard calculations.

These custom results are displayed on the monitor, logged with results, and included on hard copy print-outs as if they were part of the original application.

An example of a customer generated custom equation.

$$S = \frac{1}{\pi^2} \left| \int_0^{2\pi} \int_0^1 \exp(2\pi i \Delta W(\rho, \theta)) \rho d\rho d\theta \right|^2$$

Custom results with statistics

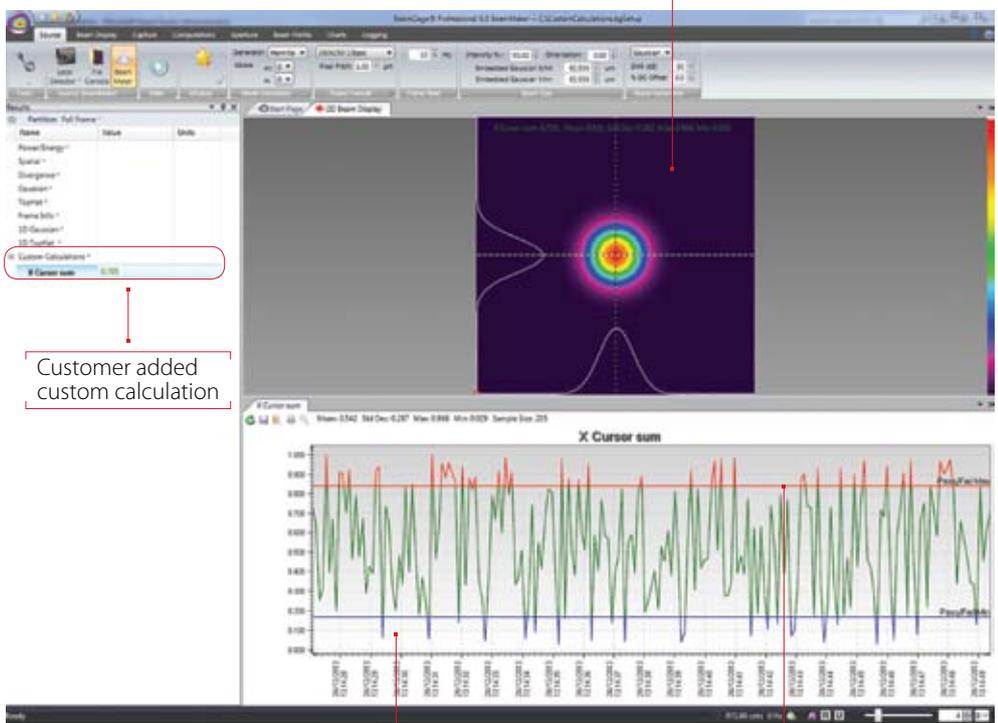

Customer added custom calculation

Custom results being plotted

Custom results with pass/fail turned on

3.2.1.3 Software Comparison Chart

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
Features Overview	<p>User selectable for either best “accuracy” or “ease of use”</p> <p>Supports our patented Ultracal algorithm plus Auto-setup and Auto-exposure capabilities</p> <p>Extensive set of ISO quantitative measurements</p> <p>Support for USB, GigE and Pyrocam III and Pyrocam IV cameras.</p> <p>New Beam Maker® beam simulator for algorithm self validation. See below for more detailed description.</p> <p>Simultaneous 2D and 3D displays</p> <p>Multi-instance, multi-camera use</p> <p>Results synchronised to select models of Ophir power/energy meters. Supported products include: Vega, Nova II, Pulsar, USBI and Juno, in both 32 and 64bit OS. (Quasar is not supported)</p> <p>Supports Satellite windows on multiple monitors</p> <p>Continuous zoom scaling in both 2D and 3D</p>	<p>Supports InGaAs and large format L11059 cameras</p> <p>Window partitioning to allow analysis of multiple beams from a single camera image</p>
	<p>Camera ROI support on USB and Firewire cameras</p> <p>Manual and Auto-aperturing to reduce background effects</p> <p>Pass/Fail on all results items, w/multiple alarm options</p> <p>Beam Pointing Stability scatter plot and stripchart results</p> <p>Full featured logging capabilities in a reloadable industry standard data file format</p> <p>Configurable Report Generator that allows cut and paste of results, images and settings.</p>	<p>NET Automation interface that allows for remote control. Examples in LabView, Excel and .Net VB</p>
Quantitative Calculations; Basic Results	<p>Supports English, German, Japanese and Chinese</p> <p>Windows OS in both 32 and 64bit.</p> <p>Multilingual GUI in English, Japanese and Chinese.</p>	
Power/Energy Results	<p>(per ISO 11145, 11146-1/-3, and 13694)</p> <p>Total power or energy (Can be calibrated or sync'd to an external Ophir power/energy meter)</p> <p>Peak power/energy density</p> <p>Min. Fluence</p> <p>Average pulse power</p> <p>Peak pulse power</p> <p>Device efficiency</p> <p>% in Aperture</p>	
Spatial Results	<p>Peak and Centroid locations</p> <p>Beam width</p> <ul style="list-style-type: none"> ■ Second Moment (D4s) ■ Knife Edge 90/10 ■ Knife Edge (User selectable level) ■ Percent of Peak (User selectable) ■ Percent of Total Energy (User selectable) ■ Encircled power smallest slit @ 95.4 ■ Moving slit (User selectable) <p>Beam diameter</p> <ul style="list-style-type: none"> ■ Average diameter (based on x/y widths) ■ Second Moment (D4s) ■ Encircled power smallest aperture 86.5 ■ Encircled power smallest aperture (User selectable level) 	

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
	Elliptical Results <ul style="list-style-type: none"> ■ Elliptical orientation ■ Ellipticity ■ Eccentricity Distance Measurement <ul style="list-style-type: none"> ■ Cursor to Crosshair ■ Centroid to Crosshair Area Results <ul style="list-style-type: none"> Beam cross-sectional area 	
Divergence	Focal Length method Far-field two-point method Far-field Wide Angle method	
Gaussian Fit	2D whole beam fits <ul style="list-style-type: none"> 1D line fits Height Width X/Y Centroid Goodness of fit Roughness of fit 	
Tophat Results	2D and 1D <ul style="list-style-type: none"> Flatness Effective Area Effective Power/Energy Fractional Effective Power/Energy Effective Average Fluence Uniformity Plateau Uniformity Edge Steepness 1D or 2D surface inclination 	
Other Quantitative Items	Frame Averaging Frame Summing Frame Reference Subtraction Image Convolution Camera signal/noise calculator Row and Column summing with results loggable	Scalable Intensity Histogram, exportable X or Y axial off axis image correction
Beam Stability Displays and Results	(per ISO 11670) <ul style="list-style-type: none"> Pointing Stability of Centroid <ul style="list-style-type: none"> ■ Scatter Plot display w/histogram ■ Mean Centroid ■ Azimuth angle of the scatter ■ Stability (M/m/S) ■ Max Radius ■ X/Y centroid/peak Strip chart plots ■ Sample/Time controlled ■ Pass/Fail limits ■ Auto scaling ■ Beam Width/Diameter Strip Charts with Results ■ X/Y M/m beam widths plots ■ Beam Diameter plot ■ Mean/Std Dev/Min/Max results displayed ■ Power/Energy Strip Charts ■ Total Power/Energy plot ■ Peak fluence plot ■ Avg Power plot ■ Elliptical Results Strip Chart ■ Elliptical orientation plot ■ Ellipticity plot ■ Eccentricity plot ■ Mean/Std Dev/Min/Max results displayed 	
Custom Calculations		User can program-in own set of calculations
Beam Profile Display Options	Utilizes advanced hardware accelerated graphics engines. All display windows can be satellited to utilize multiple display monitors. Can open one each simultaneous 2D and 3D beam display windows Common color palette for 2D and 3D displays Can open X and/or Y 1D beam slice profiles overlaid onto the 2D or 3D displays or in separate windows	

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
	Continuous software zooming in both 1D, 2D and 3D displays Pan to any detector location Continuous Z axis display magnitude scaling Multiple 128 color palettes user selectable Results items can be pasted into 2D, 3D, 1D, Pointing stability or Chart display windows.	Able to partition the camera imager into multiple regions with separate results.
1D Features	Available overlaid with 2D and 3D or in separate windows X any Y plots on separate or combined displays 1D displays with basic results and column row summing option Tophat 1D displays with Tophat results Gaussian 1D displays with Gaussian fit results 1D Profile display of the Gauss fit results on 1D, 2D and 3D displays	
2D Features	Continuously zoomable and resizable displays in satellitable window Continuous Z axis display magnitude scaling Zoomable to subpixel resolution for origin and cursor placements Pixel boundaries delineated at higher zoom magnifications Adjustable Cursors that can track peak or centroid Adjustable Crosshairs that can track peak or centroid Adjustable manual apertures Viewable Auto-aperture placement Displayed beam width marker Integrated Mouse actuated pan/zoom controls Separate 2D pan/zoom window to show current view in 2D beam display Manual or fixed origin placement	Ability to create partitions using the manual aperture controls
3D Features	3D graphics utilize solid surface construction with lighting and shading effects Integrated Mouse actuated pan/zoom/tilt/rotate controls Selectable Mesh for drawing speed vs resolution control Continuously zoomable and resizable displays in satellitable window Continuous Z axis display magnitude scaling User enabled backplanes with cursor projections	
Partitioning		Users can subdivide the imager into separate beam measurement regions. All enabled results are computed inside of each partition The manual aperture is used to define and create rectangular partition When partitioning is enabled some new results items will be enabled Centroid measurements between beams in each partition can be performed Partitioned imagers must have a single origin common to all partitions. All coordinate results are globally referenced to this single origin
Statistical Analysis	Performed on all measurement functions with on-screen display <ul style="list-style-type: none"> ■ Choices of intervals ■ Manual start/stop ■ Time from 1 second to 1000 hours ■ Frames from 2 to 99,999 Measurements reported ■ Current frame data, Mean, Standard Deviation, Minimum, Maximum of each calculation performed Controls integrated with beam stability results, scatter and strip chart plots	

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
File types	<p>Industry Standard HDF5 data and setup file format which are compatible in third party applications such as MatLab and Mathematica</p> <p>Math program and Excel compatible ASCII-csv results files</p> <p>Graphics in jpg file format</p> <p>Legacy file compatibility with LBA formats</p> <p>A user defined single file output that can contain settings, beam displays, beam profiles, charts, results, etc. in either .pdf or .xps file formats</p>	
Printing	<p>Images, reports, results, graphs, charts, statistics and setup information</p> <p>Option to print many frames in a single operation</p> <p>WYSIWYG images</p>	
Pass/Fail	<p>Set Maximum/Minimum limits on all calculations and statistics</p> <p>Red/Green font color indication on result items</p> <p>Multiple choices for indication of failed parameters, including TTL pulse for external alarm</p> <p>Master pass/fail which triggers alarm on any failure</p> <p>USB signal, beep, stop, and log alarm options</p>	
Logging	<p>Video Data Logging Formats: HDF5, ASCII-csv</p> <p>Results in ASCII-csv</p> <p>Pictures 2D and 3D in jpg, gif, tiff, bmp, png file formats</p> <p>Charts in ASCII-csv</p> <p>Cursor Data in ASCII-csv</p> <p>Row/Column summed in ASCII-csv</p> <p>Continuous Logging</p> <p>Time Interval Logging</p> <p>Frame Count Logging</p> <p>Periodic Sampling</p> <p>Pass/Fail Sampling</p> <p>Burst Sampling, after a user specified time interval, sample a user specified number of frames</p>	
Exporting	<p>Convert frame buffer data to third party format</p> <p>Export a user specified number of frames from the buffer</p> <p>Export Image Data: ASCII-cvs</p> <p>Export Results: ASCII-csv</p> <p>Export Picture: jpg, gif, tiff, bmp, png file formats supported</p> <p>Export Cursor Data: ASCII-cvs</p> <p>Export Row/Column summed: ASCII-cvs</p> <p>Export Image Data in Aperture</p>	
Automation Interface (.NET)		<p>Automation Interface with examples in LabVIEW, Excel and Net VB</p> <p>Automate launch and termination of the application</p> <p>Automate start, stop, Ultracal, Auto-X and Auto Setup</p> <p>Automate the loading of application setups</p> <p>Automate control of most camera settings</p> <p>Automate a subset of the application features and controls</p> <p>Automate the capture of Binary Video Data</p> <p>Automate the acquisition of application results</p> <p>Automate the acquisition of application Images</p>
Integrated Help	<p>PDF Operators Manual</p> <p>Context Sensitive (Whats this?) Help</p> <p>Context Sensitive Hints</p>	
Signal Conditioning for Enhanced Accuracy	<p>Spiricon's patented Ultracal enables more accurate beam measurement and display. Ultracal takes a multi-frame average of the baseline offset of each individual pixel to obtain a baseline accurate to approximately 1/8 of a digital count. This baseline offset is subtracted from each frame, pixel by pixel, to obtain a baseline correction accurate to 1/8 digital count. Spiricon's Ultracal method retains numbers</p>	

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
	less than zero that result from noise when the baseline is subtracted. Retaining fractional and negative numbers in the processed signal can increase the beam width measurement accuracy by up to 10X over conventional baseline subtraction and clip level methods. Spiricon's Ultracal conforms to the best method described in ISO 11146-3:2004.	
Frame Averaging	Up to 256 frames can be averaged for a signal-to-noise ratio, S/N, improvement of up to 16X (Noise is averaged up to 1/256th [8 fractional bits]). Data is processed and stored in a 32bit format.	
Frame Summing	Up to 256 frames can be summed to pull very weak signals out of the noise. Due to the precise nature of Ultracal baseline setting, (i.e., a retention of both positive and negative noise components) summing of frames can be performed without generating a large offset in the baseline.	
Convolution (Adjacent Pixel Averaging)	Choice of 5 convolution algorithms for spatial filtering for both display and calculations. Spatial filtering improves the visual S/N.	
Beam Maker®	Beam Maker is a new feature that allows the user to model both Laguerre-Gaussian and Hermite-Gaussian laser beams in various modal configurations. With these models you have verification and validation tools that allows not only OSI but also the end user to verify BeamGage's basic beam width measurement algorithms. It can also be used to model laser beams with special input conditions such as signal-to-noise, background offset, and bits per pixel resolution. This allows the user to better understand the accuracy of measurements made under both optimum and adverse conditions. This tool provides the user with a method to validate algorithms against current ISO standards and methods. It can also be used to validate third party algorithms by making the output data available for use in third party applications.	
Camera Features	Camera features are governed by the capabilities of the various cameras that will interface with these software products, and second by which of these camera features are implemented in the software. This section will describe typical camera features supported in the application. Black Level Control (used by Ultracal and Auto-X and Auto-setup) Gain Control (used by Auto-X and Auto-setup) Exposure Control (used by Auto-X and Auto-setup) User Programmable ROI Pixel Binning Pixel Sampling Bits per pixel setting External Trigger Input Trigger Delay Strobe Output Strobe Delay External Trigger Probe Internal Trigger Probe	
Camera related features in the applications	These are features related to but not generally dependent upon the camera design. Gamma Correction Gain Correction Bad Pixel Correction Lens Applied Option Pixel scale settings Magnification settings Frame buffer settings Ultracal Enable Auto-X (auto exposure control) Perform an Auto-Setup 8/10/12/14/16 bits per pixel Select Format or ROI Measure S/N ratio	
Trigger, Capture and Synchronization Methods	Capture methods are features related to the application while Synchronization methods relate more to the abilities of the specific camera. NOTE: Frame capture rates are determined by many factors and are not guaranteed for any specific operating configuration.	

Features	BeamGage® Standard	Upgrade to BeamGage® Professional to include: (all features in Standard plus)
	<p>Trigger modes</p> <ul style="list-style-type: none"> ■ CW - captures continuously, see Capture Options below ■ Trigger-In from laser: Trigger pulses supplied to the camera ■ Strobe-Out to laser: Strobe pulses output from the camera ■ Video Trigger: Frame captured and displayed only when the camera sees a signal greater than a user set level <p>Capture options</p> <ul style="list-style-type: none"> ■ Capture options are redefined and are approached in a different manner than older products. The items listed below will allow for all of the previous methods but with more flexibility than ever before ■ Results Priority: Results priority will slow the capture rate to be in sync with the computational results and display updates ■ Frame Priority: Frame priority will slow results and display updating to insure that frames are collected and stored in the frame buffer as fast as possible (replaces block mode) ■ Stop After: Will collect a set number of frames and then stop (replaces Single-Shot mode) ■ Periodic: Will collect frame at a programmed periodic rate ■ Periodic Burst: Will collect frames in a burst at programmed periodic rates <p>Post processing is still available but is done via a different mechanism and is limited to only data file sources</p>	
Video Playback	<p>Video playback, post processing and post analysis</p> <p>User customizable playback rates</p> <p>Video file quick pan/search controls</p> <p>Whole video file playback looping with sub-selection looping</p> <p>Playback Video produced by logging</p> <p>Almost all measurements can be performed on video files</p>	
System Requirements	<p>PC computer running Windows7 (32/64) Laptop or Desktop.</p> <p>Consult factory for latest 64bit camera availability.</p> <p>Not all cameras run in all Microsoft OS versions, see camera section for specifics</p> <p>GHz Pentium style processor, dual core recommended</p> <p>Minimum 2GB RAM (4GB required for L11059 camera)</p> <p>Accelerated Graphics Processor</p> <p>Hard drive space suitable to hold the amount of video data you expect to store (50-100 GB recommended)</p> <p>Provision for PCI-Express, FireWire, USB2 or Gigabit Ethernet input depending on camera. Interface hardware provided with camera except USB2.</p>	Minimum 3-4GB RAM

Ordering Information

Item	Description	P/N
190 - 1100nm BeamGage Standard Beam Profiler Systems (camera and software)		
BGS-USB-SP503	BeamGage Standard software, software license, 1/2" format 640x480 pixel camera with 4.5mm CCD recess. Comes with USB cable and 3 ND filters.	SP90197
BGS-USB-SP620	BeamGage Standard software, software license, 1/1.8" format 1600x1200 pixel camera with 4.5mm CCD recess. Comes with USB and cable and 3 ND filters.	SP90199
BGS-USB3-SP300	BeamGage Standard software, software license, 1/1.8" format 1928x1448 pixel camera. Comes with USB 3.0 cable and 3 ND filters.	SP90375
BGS-GigE-OSI182000	BeamGage Standard software, software license, 1/1.8" format 1600x1200 pixel Gig-E camera with 17.5mm C-mount CCD recess. Comes with Cat5e cable, power supply with ext trigger adapter, Ethernet destop and laptop PCI-Express cards and 3 ND filters.	SP90267
BGS-USB3-LT665	BeamGage Standard Edition software, software license, 1 inch format 2752x2192 pixel camera. Comes with USB 3.0 cable and 3 ND filters.	SP90377
190 - 1100nm BeamGage Professional Beam Profiler Systems (camera and software)		
BGP-USB-SP503	BeamGage Professional software, software license, 1/2" format 640x480 pixel camera with 4.5mm CCD recess. Comes with USB cable and 3 ND filters.	SP90236
BGP-USB-SP620	BeamGage Professional software, software license, 1/1.8" format 1600x1200 pixel camera with 4.5mm CCD recess. Comes with USB and cable and 3 ND filters.	SP90238
BGP-USB3-SP300	BeamGage Professional software, software license, 1/1.8" format 1928x1448 pixel camera. Comes with USB 3.0 cable and 3 ND filters.	SP90376
BGP-GigE-OSI182000	BeamGage Professional software, software license, 1/1.8" format 1600x1200 pixel Gig-E camera with 17.5mm C-mount CCD recess. Comes with Cat5e cable, power supply with ext trigger adapter, Ethernet destop and laptop PCI-Express cards and 3 ND filters.	SP90268
BGP-USB3-LT665	BeamGage Professional Edition software, software license, 1 inch format 2752x2192 pixel camera. Comes with USB 3.0 cable and 3 ND filters.	SP90378
BGP-USB-L11059	BeamGage Professional software, software license, 35mm format 4008x2672 pixel camera. Comes with universal power supply, USB cable and 3 ND filters.	SP90320
1440 - 1605nm BeamGage Standard Beam Profiler Systems (camera and software)		
BGS-USB-SP503-1550	BeamGage Standard software, software license, 1/2" format 640x480 pixel camera with 4.5mm CCD recess. Phosphor coated to 1550 nm. Comes with USB cable and 3 ND filters.	SP90198
BGS-USB-SP620-1550	BeamGage Standard software, software license, 1/1.8" format 1600x1200 pixel camera with 4.5mm CCD recess. Phosphor coated to 1550 nm. Comes with USB and cable, 3 ND filters.	SP90200
BGS-GigE-OSI182000-1550	BeamGage Standard software, software license, 1/1.8" format 1600x1200 pixel Gig-E camera with 17.5mm C-mount CCD recess. Phosphor coated 1550nm sensor. Comes with Cat5e cable, power supply with ext trigger adapter, Ethernet destop and laptop PCI-Express cards and 3 ND filters.	SP90270
BGS-USB3-LT665-1550	BeamGage Standard Edition software, software license, 1 inch format 2752x2192 pixel camera. Phosphor coated 1550nm sensor. Comes with USB 3.0 cable and 3 ND filters.	SP90384
1440 - 1605nm BeamGage Professional Beam Profiler Systems (camera and software)		
BGP-USB-SP503-1550	BeamGage Professional software, software license, 1/2" format 640x480 pixel camera with 4.5mm CCD recess. Phosphor coated to 1550 nm. Comes with USB cable and 3 ND filters.	SP90237
BGP-USB-SP620-1550	BeamGage Professional software, software license, 1/1.8" format 1600x1200 pixel camera with 4.5mm CCD recess. Phosphor coated to 1550 nm. Comes with USB and cable, 3 ND filters.	SP90239
BGP-GigE-OSI182000-1550	BeamGage Professional software, software license, 1/1.8" format 1600x1200 pixel Gig-E camera with 17.5mm C-mount CCD recess. Phosphor coated 1550nm sensor. Comes with Cat5e cable, power supply with ext trigger adapter, Ethernet destop and laptop PCI-Express cards and 3 ND filters.	SP90271
BGP-USB3-LT665-1550	BeamGage Professional Edition software, software license, 1 inch format 2752x2192 pixel camera. Phosphor coated 1550nm sensor. Comes with USB 3.0 cable and 3 ND filters.	SP90385
900 - 1700nm BeamGage Professional Beam Profiler Systems (camera and software)		
BGP-USB-XC130	BeamGage Professional software, software license, 320x256 pixel InGaAs camera with C-mount recess. .9 to 1.7um spectral band. Comes with universal power supply, USB cable, external trigger cable and 3 ND filters (consult factory for other camera options).	SP90241

Ordering Information

Item	Description	P/N
13 - 355nm & 1.06 - 3000µm BeamGage Standard and one window is included, BeamGage Professional and additional windows are optional		
PY-III-HR-C-A	Pyroelectric array detector, chopped, Grade A, one Gigabit Ethernet port and BeamGage Standard included. To complete this order, you must add an Interchangeable Window part number to accompany this system (see below).	SP90364
Windows for Pyrocam IIIHR		
PY-III-HR-W-BK7-1.064	Pyrocam III-HR window assembly, BK7, A/R coated for 1.064µm	SP90365
PY-III-HR-W-SI-1.05-2.5	Pyrocam III-HR window assembly, Si, A/R coated for 1.05 to 2.5µm	SP90366
PY-III-HR-W-SI-2.5-4	Pyrocam III-HR window assembly, Si, A/R coated for 2.5 to 4µm	SP90367
PY-III-HR-W-GE-3-5.5	Pyrocam III-HR window assembly, Ge, A/R coated for 3 to 5.5µm	SP90368
PY-III-HR-W-GE-10.6	Pyrocam III-HR window assembly, Ge, A/R coated for 10.6µm	SP90369
PY-III-HR-W-GE-8-12	Pyrocam III-HR window assembly, Ge, A/R coated for 8 to 12µm	SP90370
PY-III-HR-W-ZNSE-10.6	Pyrocam III-HR window assembly, ZnSe, A/R coated for 10.6µm	SP90371
PY-III-HR-W-ZNSE-2-5	Pyrocam III-HR window assembly, ZnSe, A/R coated for 2 to 5µm	SP90372
PY-III-HR-W-BaF2-Uncoated	Pyrocam III-HR window assembly, BaF2 uncoated for 193 to 10µm	SP90373
PY-III-HR-W-POLY-THZ	Pyrocam III-HR window assembly, LDPE, uncoated for Terahertz wavelengths	SP90374
PY-IV-C-A	Pyroelectric array detector, chopped, Grade A, one Gigabit Ethernet port and BeamGage Standard included. To complete this order, you must add an Interchangeable Window part number to accompany this system (see below).	SP90296
Windows for Pyrocam IV		
PY-IV-W-BK7-1.064	Pyrocam IV window assembly, BK7, A/R coated for 1.064µm	SP90301
PY-IV-W-SI-1.05-2.5	Pyrocam IV window assembly, Si, A/R coated for 1.05 to 2.5µm	SP90302
PY-IV-W-SI-2.5-4	Pyrocam IV window assembly, Si, A/R coated for 2.5 to 4µm	SP90303
PY-IV-W-GE-3-5.5	Pyrocam IV window assembly, Ge, A/R coated for 3 to 5.5µm	SP90304
PY-IV-W-GE-10.6	Pyrocam IV window assembly, Ge, A/R coated for 10.6µm	SP90305
PY-IV-W-GE-8-12	Pyrocam IV window assembly, Ge, A/R coated for 8 to 12µm	SP90306
PY-IV-W-ZNSE-10.6	Pyrocam IV window assembly, ZnSe, A/R coated for 10.6µm	SP90307
PY-IV-W-ZNSE-2-5	Pyrocam IV window assembly, ZnSe, A/R coated for 2 to 5µm	SP90308
PY-IV-W-ZNSE-UNCOATED	Pyrocam IV window assembly, ZnSe, uncoated	SP90336
PY-IV-W-POLY-THZ	Pyrocam IV window assembly, LDPE, uncoated for Terahertz wavelengths	SP90309
Software Upgrades		
BGS TO BGP UPGRADE	Upgrade BeamGage Standard Edition to Professional Edition. Requires a new camera key to activate.	SP90233
Camera Accessories		
USB A-B Cable	USB Cable with A to B connectors, 5 meter length.	SP90204
USB A-mini B Cable	USB Cable with A to mini-B connectors, 5 meter length.	SP90205
USB-Pass/Fail Cable	Output Pass/Fail signals when BeamGage is in out mode	SP90060
Optical Trigger for SP Cameras	Optical trigger assembly which can be mounted on camera or separately to sense laser pulses and synchronize SP cameras with pulses. Comes with a BNC cable. Stand for mounting sold separately.	SPZ17005
PCI-Express USB 3.0 Desktop Card	Optional PCI-Express, USB 3.0 card for installation to desktop PCs. Compatible with all BeamGage USB 3.0 systems.	SP90386
PCI-Express USB 3.0 Laptop Card	Optional PCI-Express, USB 3.0 card for installation to laptop PCs. Compatible with all BeamGage USB 3.0 systems.	SP90387